

Bacteria Sampling Kits have been designed for collecting and preserving pathogenic bacteria specimen from nasopharynx, wound, rectum and feces. Bacteria sampling kit contains 1 flocked swab and 1 tube with liquid Amies medium. The sampling system developed by MandeLab is highly suitable for automated & conventional processing. Customization is welcome.



Bacteria Sampling Kit Contains


Disposable Sampling Swab(s)

MS-NF7801 is a Flocked Sampling Swab with PP stick used for cell & virus specimen collection, it has been designed for flu testing and nasopharyngeal clinical diagnostics. It utilizes state of the art “spray on technology” that the flocking process by means of an electro-static charge perpendicularly attaches millions of nylon microfibers on the medical grade handle tip. The flocked swab is ideal for collecting large amount of cells and rapid elution of the specimens that instantly releases the cells into the transport medium. It has been well recognized and adopted by the diagnostic test kit manufacturers who produce reagents in molecular genetics, forensics and clinical laboratories sectors. The perpendicular nylon fibers act like a soft brush, allowing the improved collection and release of both cellular and liquid samples.

Amies Transport Medium

The liquid Amies transport medium is mainly composed of Sodium Chloride, Disodium Hydrogen Phosphate, Sodium Thioglycolate and HEPES. It is used to preserve and transport the collected bacteria samples, and only for in vitro analysis and testing, not for therapeutic use.

How to Use the Bacteria Sampling Kit?

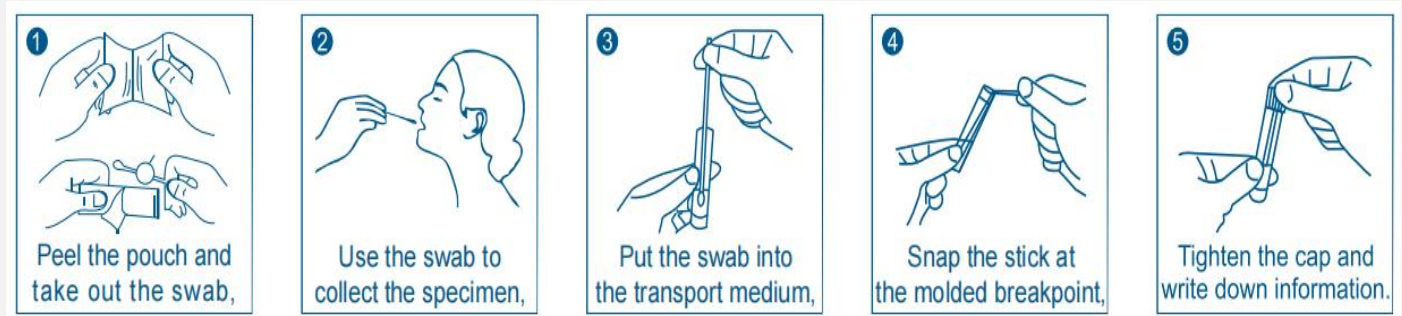

Storage and Disposal

The optimum storage temperature is 2°C to 30 °C. Do not incubate, or freeze prior to use, improper storage will result in a loss of testing. Medical waste must be disposed of in compliance with the local legislation. Please take the appropriate precautions for infected material if necessary.

Product Pictures

